

Manual

EPIC: Elution Profile-based Inference of Protein Complex Membership

Welcome to the EPIC installation manual, in this document we will have a step-by-step guide on installing, and running the EPIC pipeline.

1 Installing dependencies

In order to run EPIC locally on your machine docker needs to be installed, and for easy installation of the EPIC docker image we recommend installing kitematic, as well as creating a docker hub login. Additionally to the docker suit we also recommend installing the current version of Cytoscape.

Docker:

For mac OS X follow [this](#) instructions

Kitematic:

The current version of Kitematic can be downloaded and installed from [here](#).

Cytoscape:


Cytoscape is available from [this](#) website.

2 Building the docker image

The docker EPIC image is online available and can easily installed by running the following command line in the Terminal/shell:

```
$ docker pull baderlab/bio-epic
```

Alternatively the docker can also be installed using the Kitematic user interface as follows:

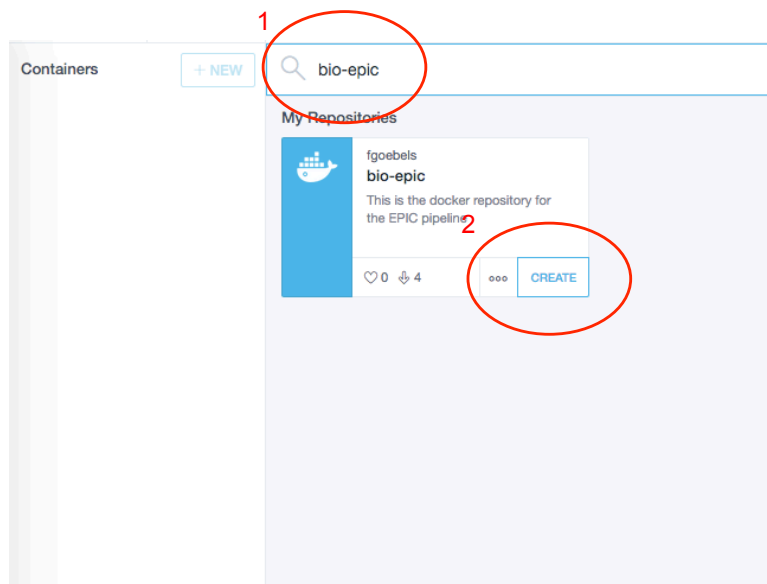

1. First starting Kitematic and entering bio-epic into the search mask
2. Clicking on the CREATE button on the Repository selection

3 Input files

EPIC accepts three types of input files:

3.1 Elution profile files

This is a tab separated file contains the elution profile for all the proteins in one distinct co-fractionation experiment, thus in case of multiple experiments each experiment has its own elution profile file. The first line of the elution profile file is the header line where the first column is ignored, and each subsequent column contains the name of the fraction. Each following line of the elution profile file contains the elution profile for one protein, where the first column contains the protein ID and the following columns each contain the MS value (either spectral counts or intensity) of the protein in each distinct fraction.

3.2 Reference protein complexes (optional)

This file can be optionally supplied as protein complex reference set, instead of having them automatically generated from STRING, IntAct, and GO. This file does not have a header line, and each row of the file describes one protein complex by concatenating all its member with tab characters.

3.3 Functional annotation network (optional)


This file can be supplied as an alternative source for functional annotation data that should be used by the EPIC pipeline. This file is also tab separated and has a header line that contains the name of each supplied functional annotation score. Each row contains the scores for each protein pair, where the first two columns contain the protein IDs, and the following the annotation scores for the protein interaction.

4 Input folder structure

For EPIC to properly run the user needs to create a folder that will be linked with the docker image, from which the docker can read and write data from the user. This folder should contain a sub-folder for each elution project the user wants to process, and those folder store all the elution profile files for the given project. In

turn EPIC will read in the data and once done will create an output folder in the same directory that contains all the output files and has the suffix “inputfoder_out”.

For example, here we created a folder named input, which contains a folder named test, which in turn contains 4 elution profile files.

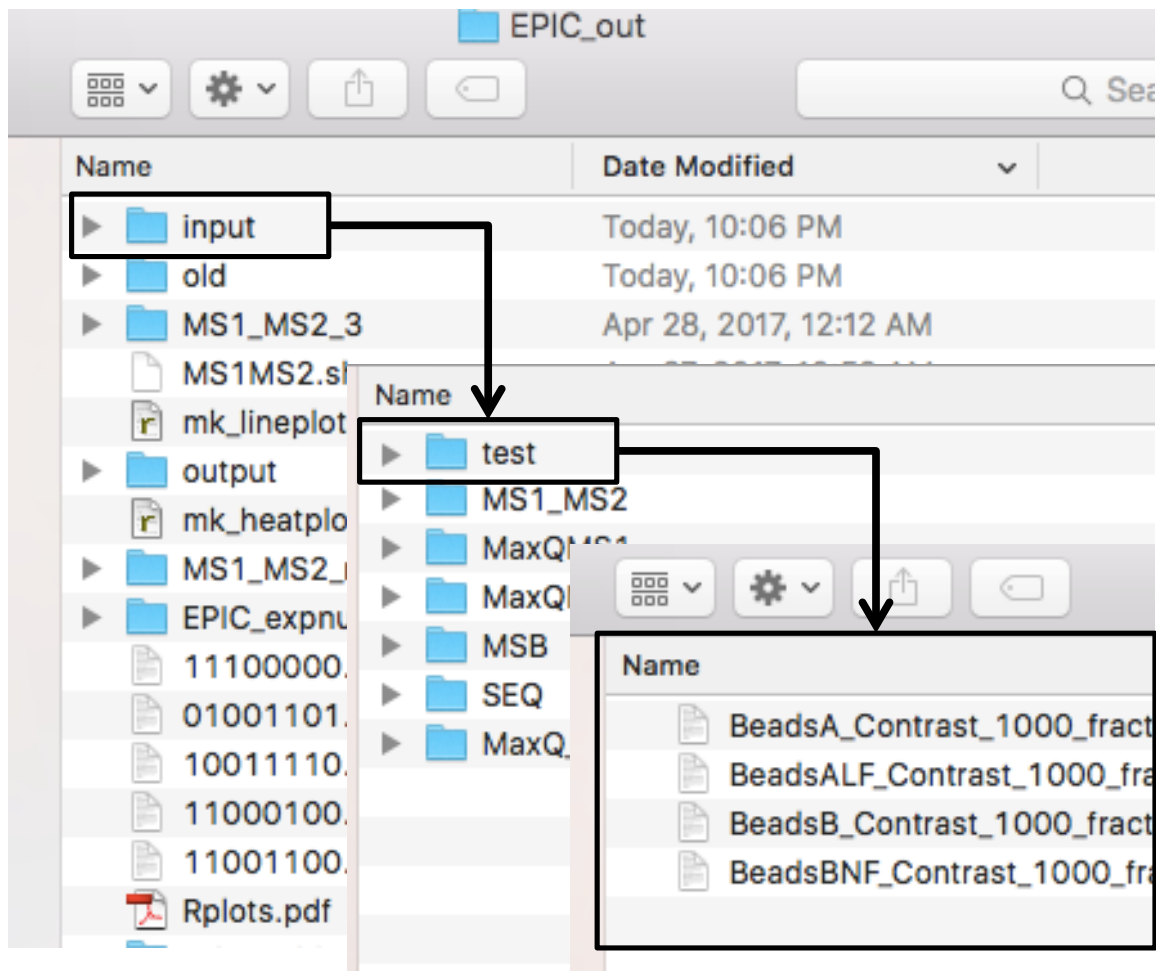

5 Starting EPIC

For mac OS X and Jupyter EPIC can be started by opening a terminal window and first using `cd` to change to the input folder (see 4) and then using the following command:

```
DIR="$( cd "$( dirname "${BASH_SOURCE[0]}" )" && pwd )"
docker run --add-host="localhost:${(ifconfig en0 | grep inet | grep -v inet6 | awk '{print $2}')} -it --rm -p
8888:8888 -v "$DIR:/home/jovyan/work/input" baderlab/bio-epic start-notebook.sh --
NotebookApp.iopub_data_rate_limit='100000000' --NotebookApp.token=""
```

After starting the EPIC container several messages will appear and one of the messages will be the Jupyter token that will allow you to connect to the Jupyter notebook.

```
[I 02:51:13.111 NotebookApp] Use Control-C to stop this server and shut down all kernels
C 02:51:13.112 NotebookApp]

Copy/paste this URL into your browser when you connect for the first time,
to login with a token:
http://localhost:8888/?token=7a566b73b38a076c8ba5dba7120214a4a9a35af09de67fd1
```

The squared red box shows the token and in the next step you need to copy paste this token into your web browser of your choice. Which will lead to this:



Form here on clicking on EPIC.ipynb will start the EPIC notebook and will result into this tab being opened:

EPIC: Elution Profile-based Inference of Protein Complex Membership

Welcome to the EPIC predictor Jupyter web service. Please use the Kitematic file organizers to upload the selected elution profiles. The following parameters need to be set in order to run EPIC.

Input directory

When starting the docker you linked one of your local folders with this docker image. This pipeline considers this folder as the input and output folder. In this folder EPIC expects a folder for each project, where each project folder contains one elution profile file for each co-fractionation experiment. The elution profile file is a tab separated file where the first column contains a protein ID and each subsequent column contains protein value for each fraction, and the file has one row for each protein. Additionally the file has an header line which contains the name of each fraction.


Feature selection:

Please select which co-elution features should be used to generate the co-elution network. We recommend using MI, Bayes, PCCN, and Apex.

6 EPIC output

After following the prompts of the Jupyter notebook several output files and plots will be generated as well as a Cytoscape network will be generated in a running Cytoscape instance.


6.1 Classifier performance plots

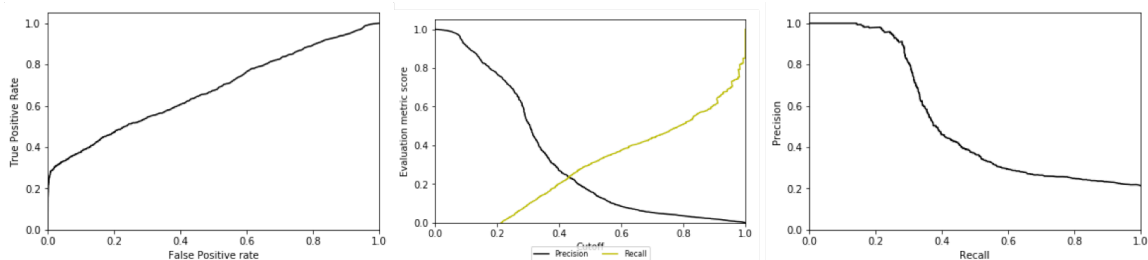

From left to right: ROC curve, Precision/Recall vs confidence score, and PR curve.

The first three plots that are generated show the classifiers performance in predicting co-complex membership from the reference data. The first plot show both precision and recall for various classifier confidence values, and this plot will give the user insight to evaluate how confident an predicted interaction with a certain score is. The following two plots show the ROC and PR curve for the classifier.


6.2 Cytoscape visualization within the browser

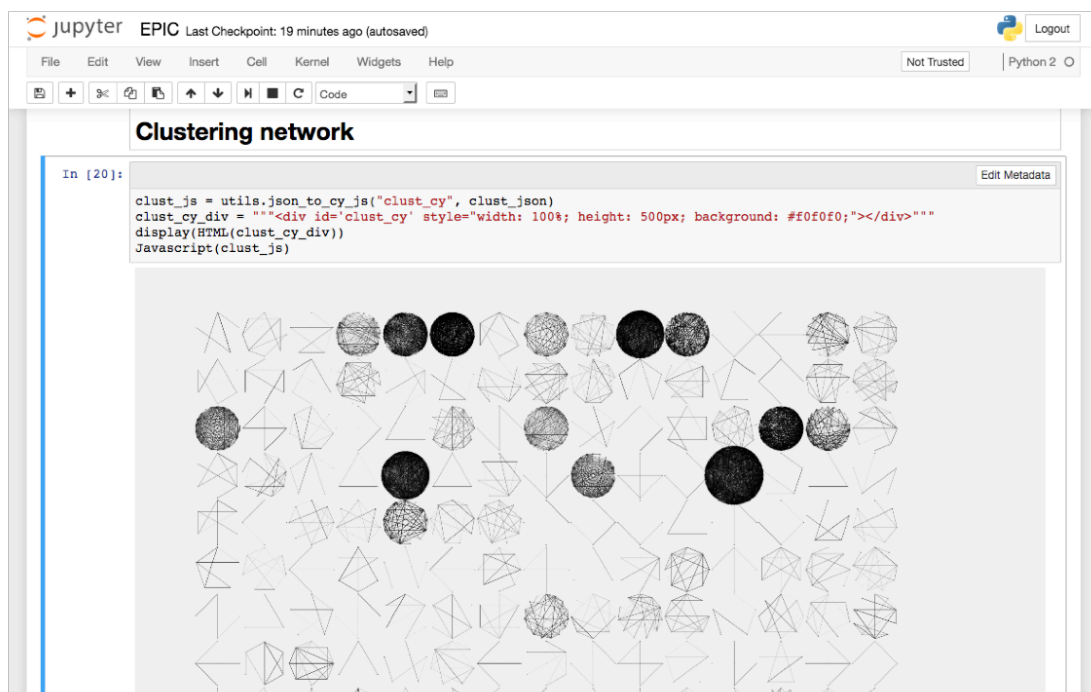

The next figure shows the interactive cytoscape.js widget that visualizes the generated protein complexes, and the user can zoom in and out as well move clusters in side this widget. This allows users to have an initial glance at their created protein complexes, and can always be generated even if the user does not have Cytoscape locally installed.


6.3 Cytoscape

Running the last cell in the EPIC Jupyter script will send the created clusters to a running Cytoscape instance and generate a new network named “EPIC clusters”, and once selecting the running Cytoscape session the user is prompted to create the view. Once the view has been created the user will see a network that contains a connected component for each predicted cluster.

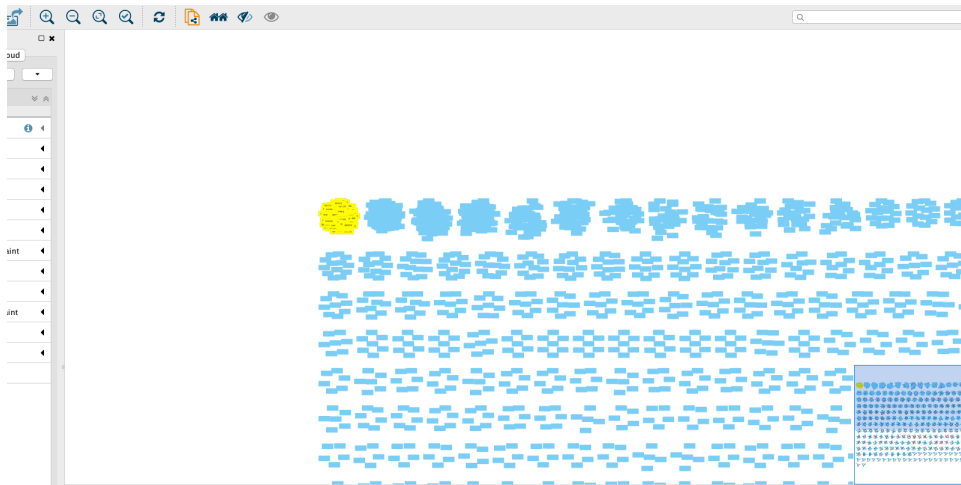


If the user has installed the clustermaker app from the Cytoscape app store he can visualize the elution profile for any selected protein as follows:


First select the JTree heatMapView from Apps:

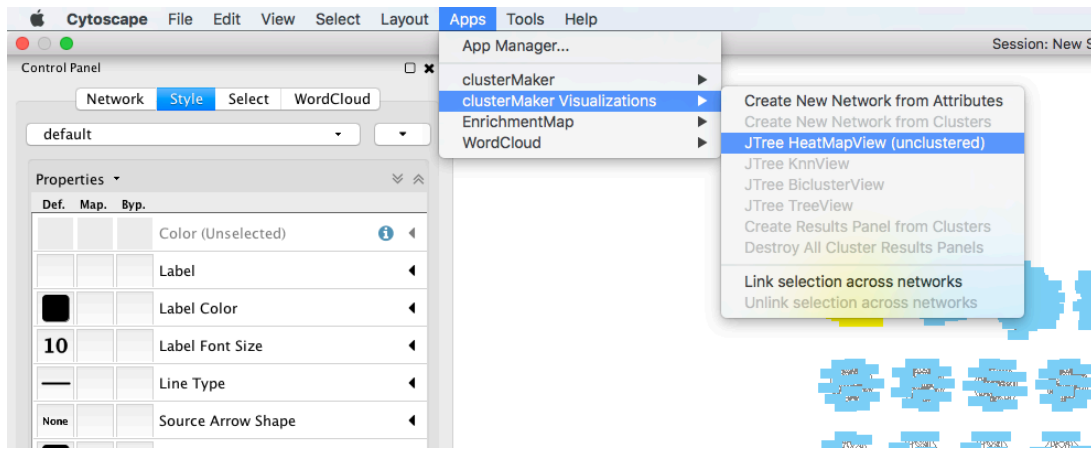

This will open the following window in which the user should select all elution score Node attributes, and check the lower box that says “Use only selected nodes/edges for clustering”:

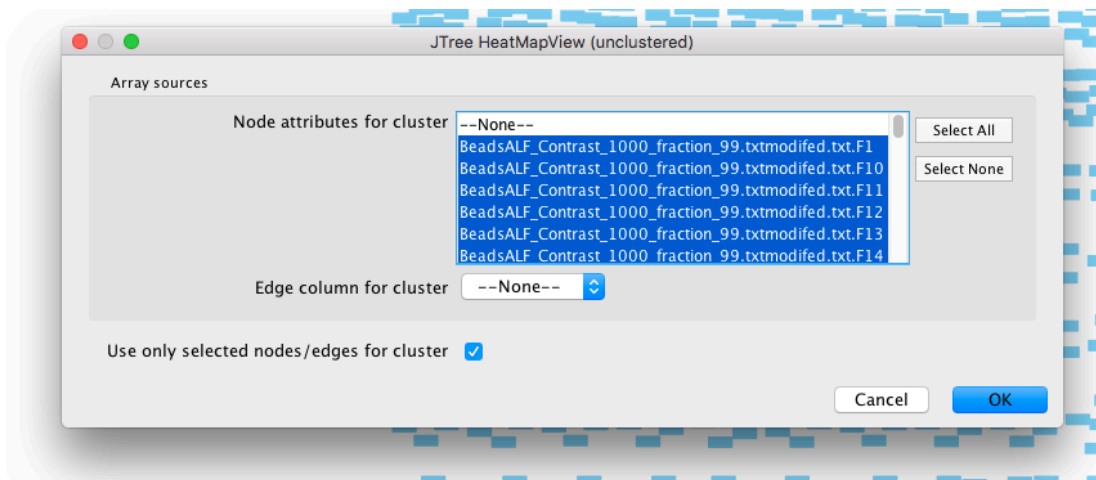

This will result in opening another window in which the user can freely browse the elution profiles of the selected proteins:

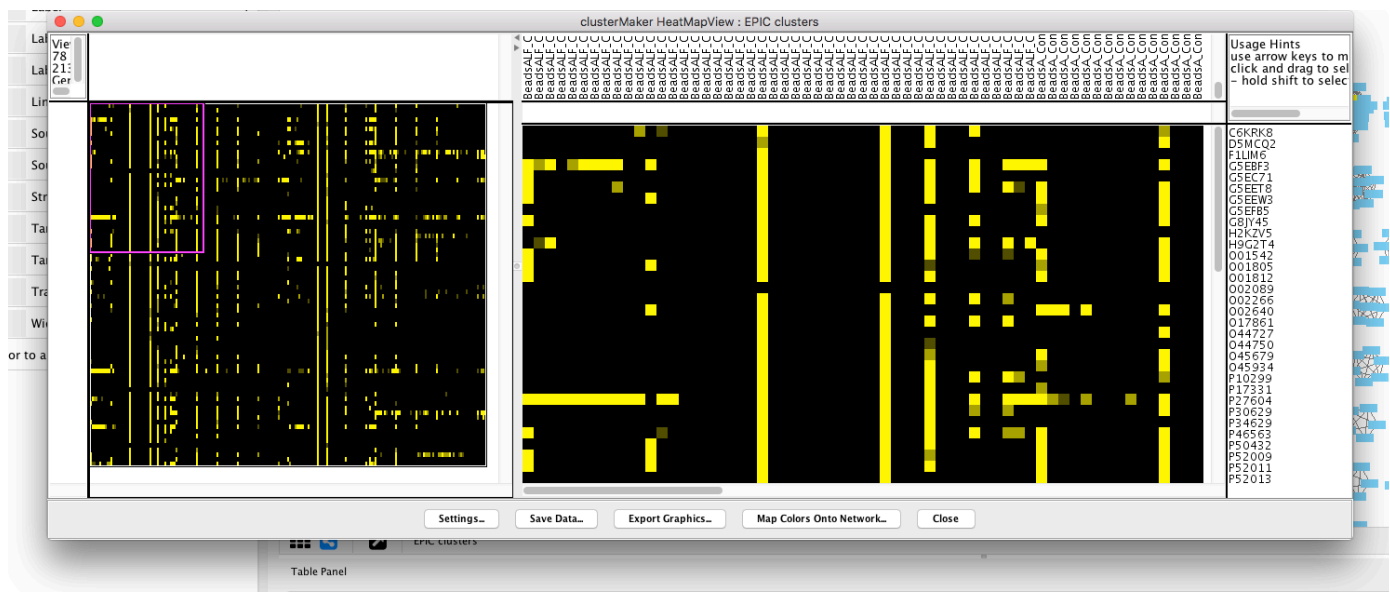

7 Collecting the output

EPIC will create an output directory in the shared folder that contains the project folder and stores several output files in it. Each file has a prefix consisting of the used classifier method and the mode in which EPIC was run with (modes are EXP, FA, COMB, and BR). For example if EPIC was run with a random forest classifier and experiment only each file will have the prefix “Out.rf.exp”. The most important files are as follows:

***.scores.txt**

contains the raw co-elution score.

***.pred.txt**

This tab-separated file contains the generated binary protein interaction network, where the first two columns stores the protein ids and the third column is the classifier confidence score.

***.comp.txt**

This is the predicted protein complex file, where each line describes one protein complex.